


10-Y73-TL or 35-Y73-TL for either 25mm or 35mm *SmartShutter* installation. Alternatively, a TLED-Y can be used in this location as an LED light source

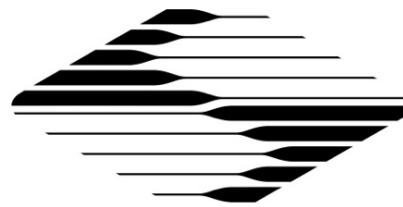

MD platforms are available for both the IX83 and IX73:
MD-83-1 = 1 sided IX-83
MD-83-2 = 2 sided IX-83
MD-73-1 = 1 sided IX-73
MD-73-2 = 2 sided IX-73

LG-Y73 collimating adapter replaces the RFA epi-illumination arm.

If the RFA is desired:
The LG-Y73-LF collimating adapter can mount to the end of the RFA. The 10-Y73-EC and 32-Y73-EC will allow connection of the *SmartShutter* or excitation wheel to the microscope directly

LB10-W32-Y73 emission wheel can be installed in one of the deck position. This wheel uses 32mm emission filters

10-Y73-EM can be used to mount a LB10-NWE emission filter wheel

